## Neck Massager NM20

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# Content

---

Important information	16
Specifications	19
Device and controls	20
Use	20
Neck Massaging	21
Massaging other parts of the body	21
Cleaning	22
Warranty	23

## Important information

---

- Only use the power supply unit supplied for the massager.
- Keep the power supply unit, cable and appliance well away from heat, hot surfaces, moisture and liquids.
- Do not touch the power supply unit when you are standing in water and always take hold of the plug with dry hands.
- Do not take hold of a device if it falls into water. Pull out the power plug and the power supply unit from the socket immediately.
- Connect the device only in that way, so that the power plug is freely accessible.
- Always switch off the appliance after use with the buttons on the control unit and unplug the power supply unit.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the device on the cable or the power supply unit.
- You must not continue to use the unit if the lead or the mains adaptor are damaged. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the unit in to the service centre.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must be supervised to ensure that they do not play with the device.

- Do not use this unit to supplement or replace medical treatment.
- Chronic complaints and symptoms may be made worse.
- Do not use the massager, or only after consulting your doctor, if:
  - you are pregnant
  - you have a cardiac pacemaker, artificial joints or electronic implants.
  - you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- Do not use the appliance near your eyes or near any other sensitive areas.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Take care when using the heating function. The appliance has hot surfaces. Do not use the appliance on small children, vulnerable individuals or people who are not sensitive to heat.
- Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

**The device is not to be used for commercial use or medical purposes.**

**If you have health concerns, consult your doctor before using the device.**

- The device must only be used for its intended purpose as described in the instruction manual.
- Using the unit for any other purpose invalidates the warranty.
- Only use the appliance in the proper place, as described in this instruction manual.

- Do not use the appliance while you are asleep or when you are lying in bed.
- Do not use the appliance before going to sleep. The massage has a stimulating effect.
- Never use the appliance while driving a vehicle or operating a machine.
- The maximum operating time for one session is 15 minutes. Prolonged use reduces the service life of the appliance.
- Very long periods of use can lead to overheating. Always allow the appliance to cool down before using it again.
- Never put the unit down next to an electric stove or other sources of heat.
- Never leave the unit unattended if it is connected to the power supply.
- Do not allow the massage cushion to come into contact with pointed or sharp objects.
- Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.
- Make sure, that no fingers or other parts of the body come between the rotating massage heads. **Risk of injuries!**
- Do not use an external time switch or a separate remote control system to operate the appliance.
- Before use
  - Check the power supply unit, cable and massager carefully for damage before each use. Do not attempt to use a defective unit.
  - Do not use the appliance if it, or the power cord show signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.

# Specifications

---

- Two direction
- Heating function
- Automatic standby after 15 minutes of use

Voltage	DC12V
Power	24 W
Frequency	50Hz

## Power Supply

This devices has an external power supply.

## Content of packaging


- 1 Flow Fitness Neck Massager NM20
- 1 Power supply unit
- 1 Instruction Manual

# Use

---

Thank you for purchasing the Flow Fitness Neck Massager NM20. With the Ab Shaper you have purchased a high quality product. We advise you to carefully read this instruction manual before using this product so that you can use it safely and benefit from it maximally.

## Device and controls
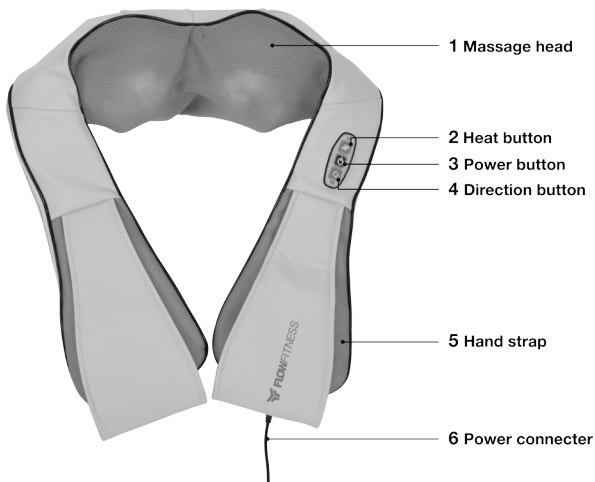
## Neck Massaging

1. Insert the power supply unit into the power connector (6)
2. Plug in the power supply unit into the socket.
3. Place the massager around the neck so that it has sufficient hold on the shoulders.
4. Guide your hands through the hand straps (5) on right and left side. By pulling on the straps you may stabilize the device in its position and adjust the massage intensity (pulling the straps stronger increases the massage intensity).
5. Press the power button (3) to start the massage.
6. Press the direction button (4) to change the direction of the massaging heads.
7. Press the heat button (2) for warming you muscles during a massage. The light in the massaging heads will indicate that the heat has been turn on. (The heating function only works with active massage function and generates a comfortable, smooth heath.) Press the heat button (2) again to turn the heat off.
8. Pressing the power button (3) again to switch the neck massager off. After 15 minutes of constant massaging the device will switch off automatically.

## Massaging other parts of the body

The device can alternatively also be used on other parts of the body, as shown in the illustration below. To apply the massager there also place your hands in the hand straps to apply pressure. You may also use the additional warming function.

Neck/Shoulder


Upper Back

Lower Back

Buttocks

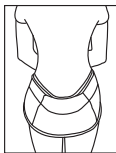Thighs

## Cleaning

---

The unit is maintenance free. If a fault does occur, just check whether the mains power adapter is properly connected.

- Before cleaning the appliance, make sure that the appliance is switched off and the power supply unit is disconnected.
- Clean the massager with a slightly damp sponge only. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol.
- Rub the massager dry with a clean dry cloth.
- Never immerse the appliance into water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

## Warranty

---

The warranty is provided by the Flow Fitness distributor of your country and the reseller where you have purchased your product. Please check <https://www.flowfitness.com> for contact information.

Neck massager NM20

[www.flowfitness.com](http://www.flowfitness.com)

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Copyright © Flow Fitness
All rights reserved. Nothing in this instruction manual may be used or reproduced in any form or way without the express written permission of Flow Fitness.